LED Light Teeth Whitening Booster

Art no: 650298

Boosts performance of teeth whitening

Claims:

- Compatible with all Beconfident® Teeth Whitening systems
- Improves up to +1 shade whiter teeth
- Works best with Beconfident® Pro Kit that contains Ferrous Gluconate
- With built in 15 min timer
- Can be used with Beconfident® original Mouth trays
- Battery included – lasts up to 20 sessions

Product description:

The Beconfident® LED Light Teeth Whitening Booster is a portable lamp using blue LED light to improve the effect of teeth whitening. It's compatible with any of the Beconfident Teeth whitening Kits, however is especially designed to work with the gel of the Pro Kit. This gel contains Ferrous Gluconate which creates a reaction for teeth whitening. Simply said, instead of generating 1 oxygen free radicals per molecules of peroxide, it generate two hydroxyl radicals, which do the same thing as the oxygen radical, except it end up with twice as many free radicals.


Active ingredients / How it works in short:

Ferrous gluconate is necessary to cause Fenton's Reaction as explained below;

In order to achieve maximum effectiveness and patient comfort, the gel has been reformulated to include ferrous gluconate for more effective peroxide activation. Ferrous gluconate provides the Iron (Fe2+) for Fenton's Reaction:


FENTON REACTION

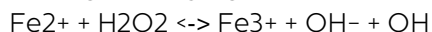In the presence of a suitable chromophoric substrate, such as the stains found on teeth, the hydroxyl radicals generated are capable of whitening teeth via oxidation of the stain's carbon molecular bonds.

Irradiating the bleaching gel with LED UV light energy even further enhances the catalytic whitening effect of Fe2+, so that a photo-Fenton process begins. The advantage of photo-Fenton chemistry is that provides enhanced amounts of bleaching free radicals via a remarkable "never-ending" regeneration reaction as shown by the reaction below:

PHOTO-FENTON REGENERATION OF FENTON'S REACTION

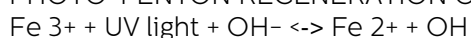Group Beconfident Europe AB

Tennisvägen 31, SE-269 33 Båstad, Sweden

E-mail: info@beconfident.info Phone: +46-(0)702 56 50 00 VAT no: SE 556853-0041-01

Corporate site: www.beconfident.info, www.beconfident.co.uk, www.beconfident.se & many more

Beconfident- Top Quality Teeth whitening and Skincare systems are designed in Sweden & Made in USA and Sweden.

Beconfident® + WhiteAmin® + SimpleSmile® are registered trademarks in Europe and Beconfident® is a registered trademark in United States of America. IOS and Android Intelligent Apps belong to Beconfident Skincare AB. Other brands are pending.

How to use:

EN: Open up the bags and take out the LED light and mouth tray. Mount the Mouth tray onto the front of the lamp. Open the battery lid on the bottom and remove the plastic film in order for the battery to become activated. To start LED light press button. Light will automatically go out after 15 minutes.

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